



DETVIEW
POWERED BY DETCOM

Admin Manual

DetCom Systems Ltd.

Version: 3.0

© 2020

www.detview.com

Content

Content.....	2
Introduction.....	2
Login page, Login.....	2
Doctors Page.....	4
Patients Page.....	5
Wrong emails page.....	8
Statistics page.....	9

Introduction


The document contains a detailed description of the admin level of the DetView system.

Login page, Login

There are additional language options available on the site, which can be selected from the drop-down menu.


Detailed description of the login:

1, You receive a message with accept link. If you clicking on this link, you agree to store your email address in the DetView system.

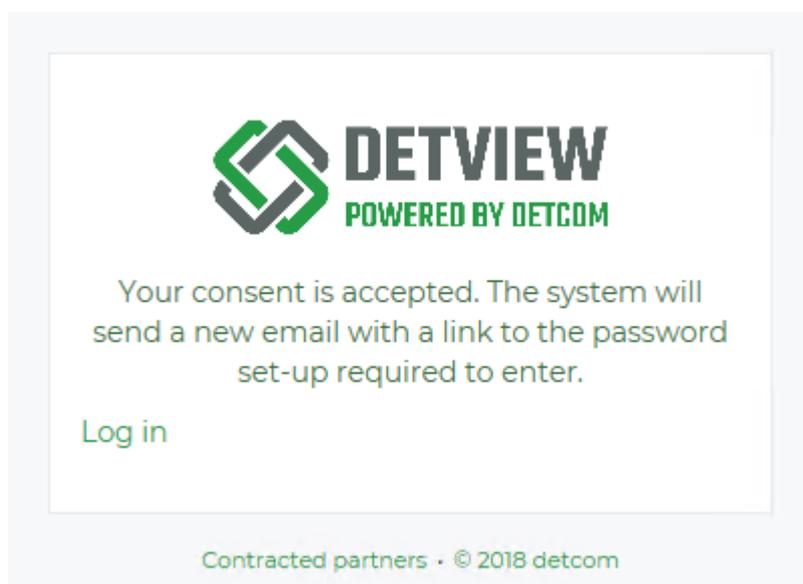


2, In the case of successful, you receive a new message with password setting link. Where you can enter own password. When you choose your password you should consider security, so it is advisable

to choose a password in the following way: minimum 8 characters, including special characters, numbers, and capitalization.

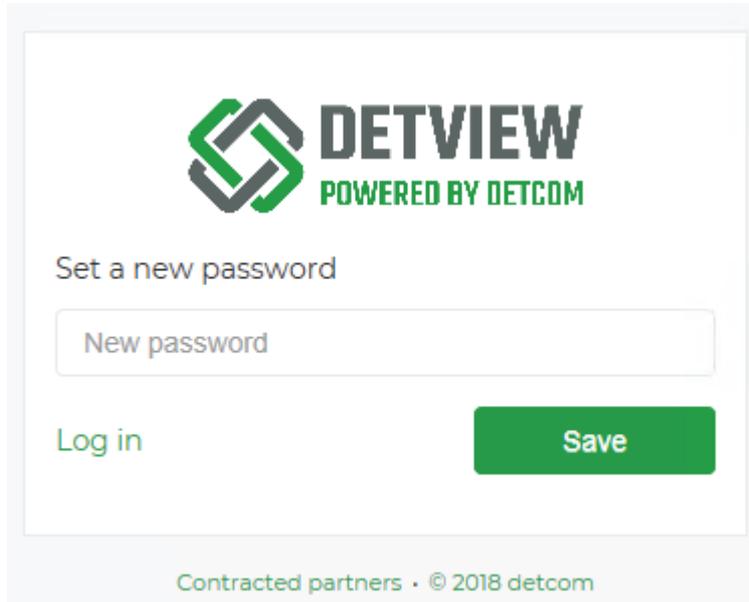

The screenshot shows a web form titled "Set a new password". At the top is the DETVIEW logo with the tagline "POWERED BY DETCOM". Below the logo, the text "Set a new password" is displayed. There is a single text input field labeled "New password". Below the input field, there is a "Log in" link and a green "Save" button. At the bottom of the form, there is a footer that reads "Contracted partners · © 2018 detcom".

3, In the case of successful save, you receive email.

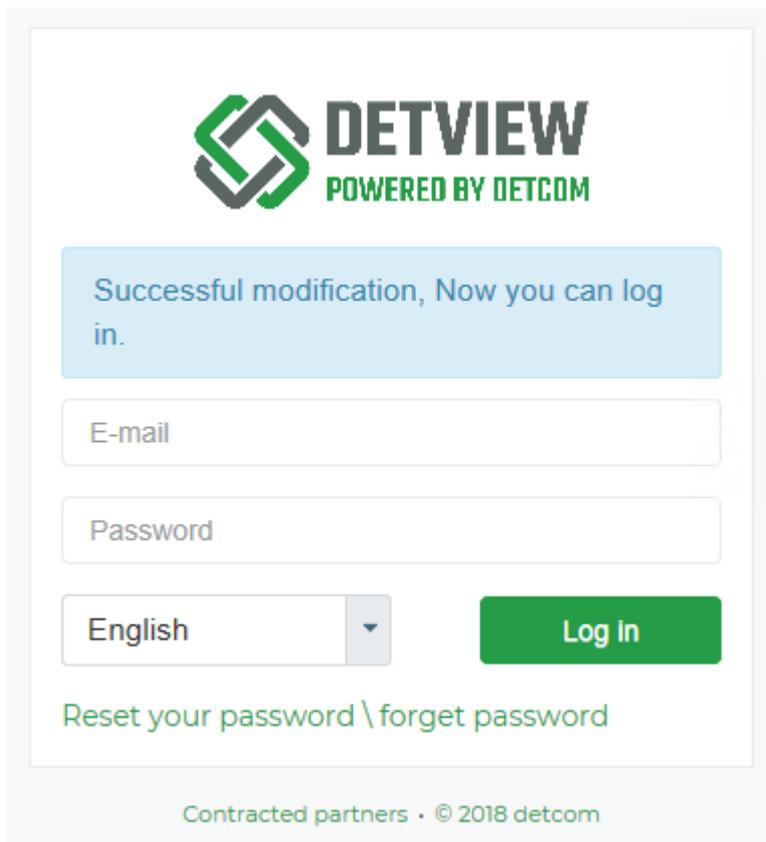

The screenshot shows a web form titled "Successful modification, Now you can log in." The message is displayed in a light blue box. Below the message, there are two text input fields labeled "E-mail" and "Password". There is also a dropdown menu currently set to "English" and a green "Log In" button. Below the input fields, there is a link that reads "Reset your password \ forget password". At the bottom of the form, there is a footer that reads "Contracted partners · © 2018 detcom".


4, Then you get your main login surface. Enter your e-mail address and password here and click on "Sign In"

Doctors Page

Here you can add or edit doctors' email address.

From the top menu, you can access the User manual, Sharing log and FAQ menus.

The system allows you to add more doctors to an email address. This may be necessary if a doctor is registered under several names in the source system.

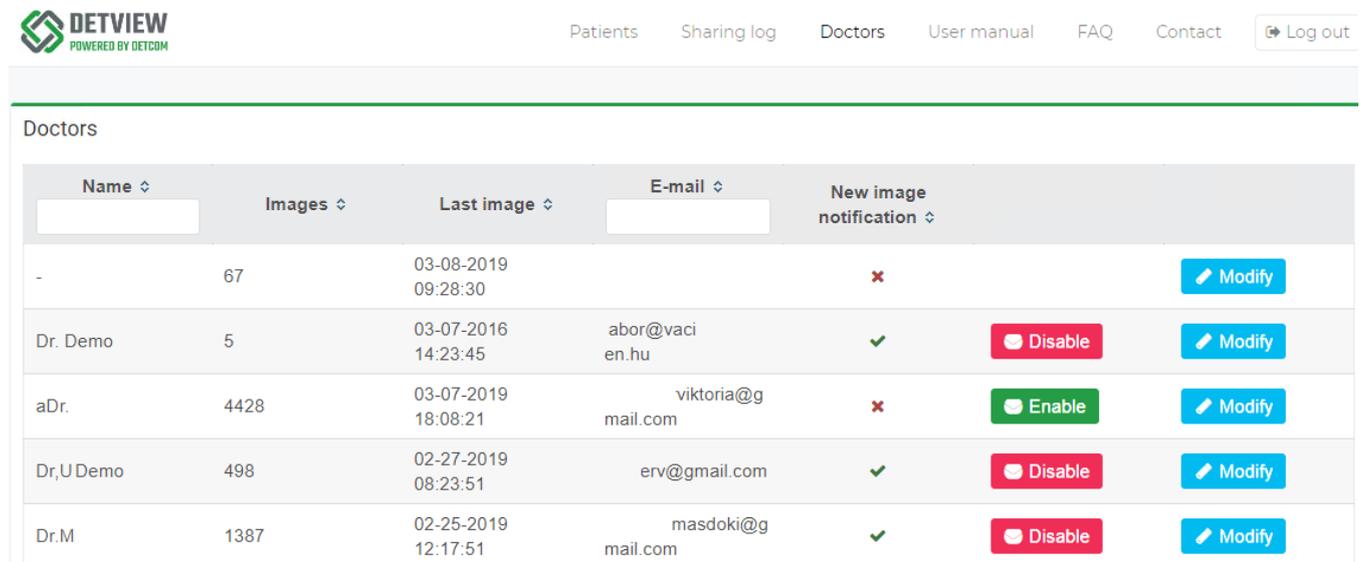

Name	Images	Last image	E-mail	New image notification	
-	67	03-08-2019 09:28:30		✗	Modify
Dr. Demo	5	03-07-2016 14:23:45	abor@vacien.hu	✓	Disable Modify
aDr.	4428	03-07-2019 18:08:21	viktorias@viktoriasmail.com	✗	Enable Modify
Dr,U Demo	498	02-27-2019 08:23:51	erv@gmail.com	✓	Disable Modify
Dr.M	1387	02-25-2019 12:17:51	masdoki@gmail.com	✓	Disable Modify

Detailed description of the page:

- Columns on the surface can be sorted as desired.
- Filling out the "Name" and the "Email" columns can result in faster searches.
- By default, only 10 doctors can be seen on the page but it can be increased to 50 by using the drop-down menu below. If you would like to see more, you can turn a page on the bottom of the page.
- By clicking on "Sign out", which is on the right, you can log out of the system.
- By clicking on the "Modify" button, you can add an email address to a particular doctor or edit the existing one. This email address determines who can access the system.

Name	Images	Last image	E-mail	New image notification	
-	0		varga.ferenc19@gmail.com	✓	Disable Modify
Dr. Demo	28	03-21-2003 07:41:28		✗	Enable Modify
dr. fer	96	10-04-2018 19:35:16	dr. fer: varga.ferenc19@gmail.com	✗	Enable Modify
doki1	1	08-29-2016 10:12:19		✗	Enable Modify
doki2	1	12-08-2016 14:39:34	doki7@mailnesia.com	✗	Enable Modify
doki3	1	06-11-2014 18:19:38	doki3@mailnesia.com	✗	Enable Modify

- If you click on the "Save" button, you save the email address and if checked in the "Send Email Notification" check box, at that moment the doctor will also receive a notification via email.
- If you click on the "Enable" button, then doctor will receive an email notification from each patient you are taking. In addition if you click "Disable" button then system will not receive email a notification email to the doctor.
- In the sharing log you can view previously shared patients the sharing doctor with assigned doctor or patients and their dates with sharing states. The green tick indicates the active shares and the red dustbin deleted shares.

Patient shares

Shared patient name	Sharing doctor	Assigned doctor	Sharing date	Date of deletion of share	Sharing status
Balázs Orbán	ferenc1980@gmail.com	92@gmail.com	09-11-2018 21:22:36		✓
Anonymous1_DICOMR M SIDEXIS	ferenc1980@gmail.com	teszt11@mailnesia.com	09-12-2018 18:39:30	09-12-2018 18:40:41	✗
Anonymous1_DICOMR M SIDEXIS	ferenc1980@gmail.com	teszt@teszt.com	09-13-2018 12:21:11	09-14-2018 07:02:23	✗
Anonymous1_DICOMR M SIDEXIS	ferenc1980@gmail.com	teszt2@teszt2.com	09-13-2018 12:47:00	09-13-2018 12:47:28	✗

Patients Page

On this page, you can download 2D and 3D images of patients assigned to the user. Here you also have the opportunity to share the patient's images with other doctors and patients.

From the top menu, you can access the User manual and if you need add-ons for the application, you can download them from here.

DETVIEW POWERED BY DETCOM						Patients	Additional download	User manual	Log out
Patients									
Doctor's office	Name	Birthday	Last image	Images					
NFRC	Alex	-15-1997	08-29-2016 08:12:19	1		1			
NFRC	Abuczki	12-29-	12-08-2016 13:39:34	1		2			
1-2 / 2 1 10									
Shared patients									
Doctor's office	Name	Birthday	Last image	Images					
NFRC	Csabáné	12-29-	12-08-2016 13:39:34	2					
1-1 / 1 1 10									

Detailed description of the page:


There are two sections on the page, one with the patients and one with the shared patients. The latter show the patients who are shared with other doctors or patients

- Columns on the surface can be sorted as desired.
- Filling the "Name" and the "Date of Birth" columns can result in faster searches.
- By clicking on the "Open" button, we get the patient's all images that can be retrieved from the system.

03-09-2020 09:05:09	Panorama Recording			
03-09-2020 10:22:29	Panorama Recording			
03-10-2020 08:07:36	3D Recording			
161-170 / 194 11 12 13 14 15 16 17 18 19 20 10				

There are several ways to download recordings on this page:

1. Clicking on the "Download image" button will immediately download the image in "jpg" format.
2. Clicking on the "Edit" button will display the recording, which can be armed as desired before downloading.

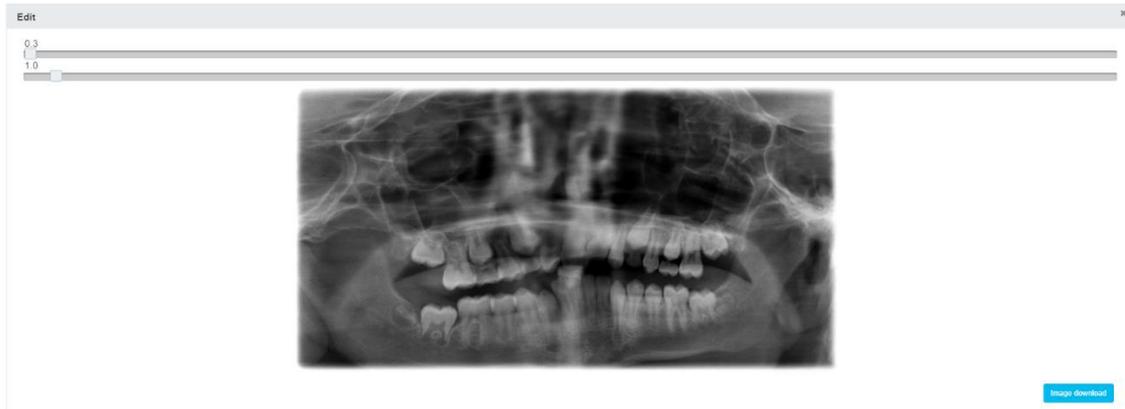


3. After clicking on the "Dicom file download" button, the recording will be downloaded in "dcm" DICOM format, so we can load it into our own SIDEXIS or other X-ray viewing systems.
 4. After clicking on the "3D without Viewer" button, the images are downloaded in DICOM "folder" format, which can also be loaded into other compatible X-ray viewing systems
 5. Under the "Download 3D with Sidexis4 Viewer" button, you can download the format that can also be written to a DVD, which also includes the X-ray viewer. (This file is the largest, more than 1GB, so the download is longer!) • If you download recordings from an older system that is not fully DICOM compatible, you will find them under the exported recordings.
- If patients has some 3D or DICOM images then you can see exported images table. Depending on the export:
 - 2D and 3D images "dcm" formats
 - 3D DICOM formats without Viewer
 - 3D DICOM formats with Sidexis4 Viewer
 - 3D formats with Sidexis3 Viewer
 - inappropriate exports are indicated by the system

Exported images		
Time ↕	Image type ↕	
11-18-2018 13:25:21	3D image with Sidexis4 Viewer	Download
10-17-2018 21:38:57	Not valid image. Please contact X-Ray center	Download
10-17-2018 21:39:56	3D image without Sidexis Viewer (DICOM)	Download
10-26-2018 08:41:57	3D image with Sidexis3 Viewer	Download
10-26-2018 08:46:32	2D image in DICOM format	Download
10-26-2018 08:48:12	2D image in JPG format	Download

1-6 / 6 << < 1 > >> 10 ▾

- If you click the "Download" button then system will make self-extracting file. Depending on the speed of the server, you have to wait for the file to be created. As soon as the executable is ready the download will start. After downloading, downloaded "exe" file name includes patient details. If you click to "exe" file then you extract it to the folder where we downloaded it. You can then launch the startup file in the unpacked file as if you were receiving a CD.
- By clicking on the "Share" button, you can share the images of the patient with a doctor or patient.

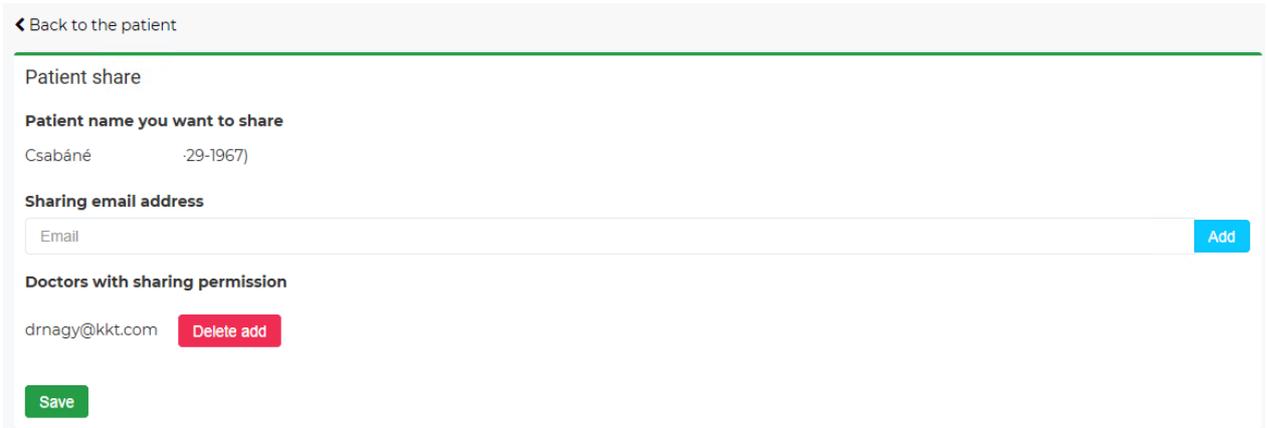

- Here enter the email address of the doctor or patient with whom we want to share the images.


First, select from the drop-down menu to share with a Patient or Doctor and enter the email address. After clicking on the "Add" button, a pop-up window will appear, where you will be informed about the data protection contribution. If you are impressed, click the "Save" button.

However, we also have the option for deleting added email addresses by clicking on the "Delete add" then "Save" button. In both cases, the system sends a notification notice to the given Email addresses.

- By default, only 10 patients and 10 images can be seen on the pages but it can be increased to 50 by using the drop-down menu below. If you would like to see more, you can turn a page on the bottom of the page.
- By clicking on "Sign out", which is on the right, you can log out of the system.

Wrong emails page

Here you can see the email addresses that the system was unable to send because they are not real email addresses or have simply been recorded incorrectly. This is a notification interface to correct the error.

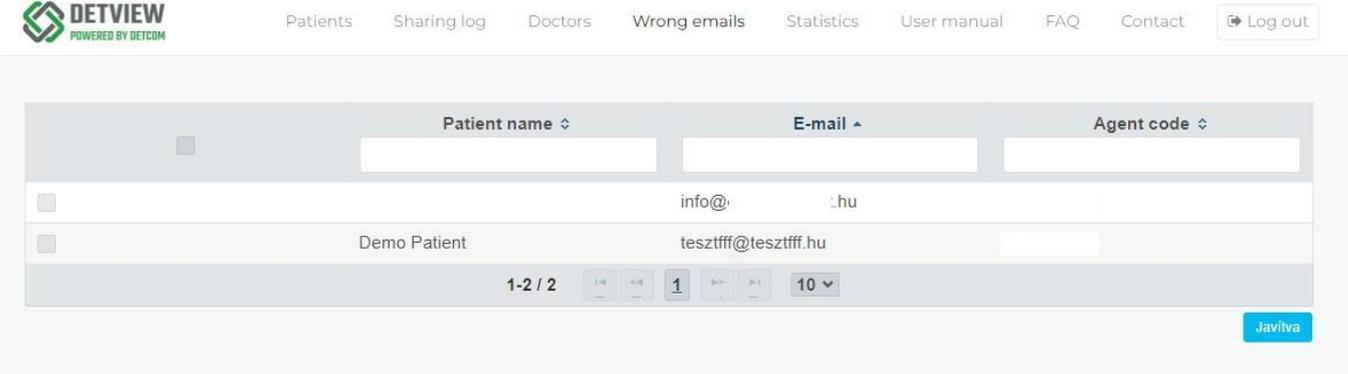

Patients Sharing log Doctors Wrong emails Statistics User manual FAQ Contact Log out

Patient name E-mail Agent code


<input type="checkbox"/>	info@	.hu
<input type="checkbox"/>	Demo Patient	tesztfff@tesztfff.hu

1-2 / 2 1 10

Javitva

Statistics page

On this page you can view your entries and downloads. Who downloaded what, when and how many times.

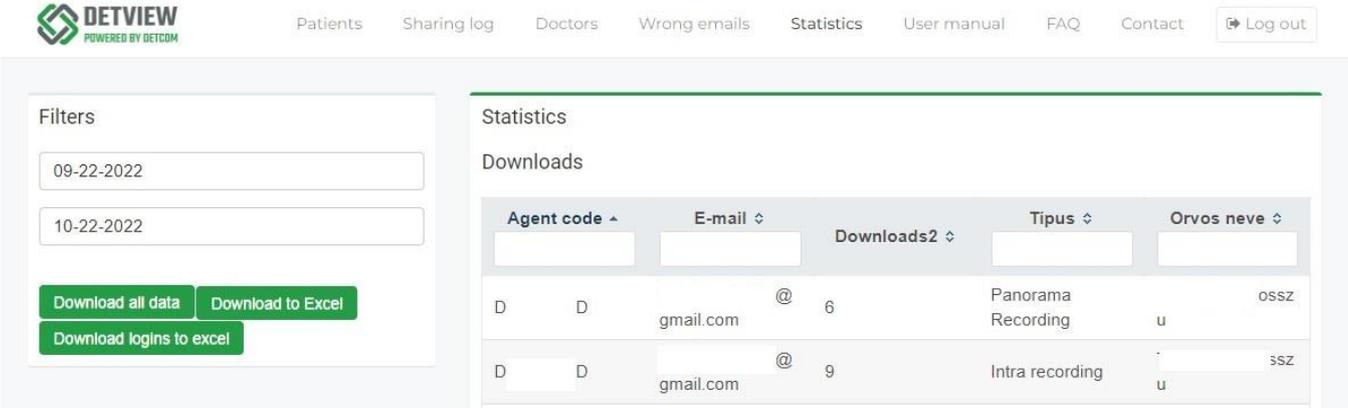

Patients Sharing log Doctors Wrong emails Statistics User manual FAQ Contact Log out

Filters

09-22-2022

10-22-2022

Download all data Download to Excel
 Download logins to excel

Statistics

Downloads

Agent code	E-mail	Downloads2	Tipus	Orvos neve
D D	gmail.com @	6	Panorama Recording	u ossz
D D	gmail.com @	9	Intra recording	u ssz

The whole database can be downloaded in excel format, so you can then set up different queries, even with graphs.

